


Buttonfix[®]

Adjustable Loop Device

One of the strongest adjustable suspensory fixation devices for cruciate reconstruction




Buttonfix[®] Adjustable Loop Device

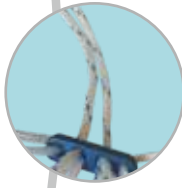


Adjustable Suture

- Enables calibrating the loop to its optimum size.
- Allows the surgeon to maximize the amount of tendon in the femoral tunnel, thereby optimizing the graft-tunnel healing interface and the chance of a successful outcome.

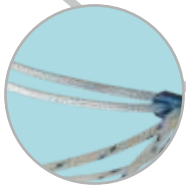


Flipping BioFiber[®] Suture

- To flip the button over the lateral cortex of femur.

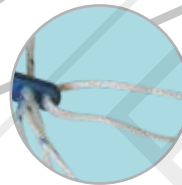


Pulling BioFiber[®] Suture

- To pull the graft from tibial and femoral tunnel.

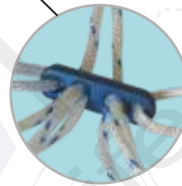

Titanium Button

- Strong with low profile.
- Compatible with common extension devices.

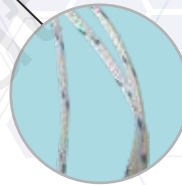

UHMWPE Suture Loop

- Provide high strength to the loop.
- As strong as Closed loop.



Graft Support Frame

- A larger and thicker surface contact with the graft.
- Protects the graft from damage and any abrasion.

Ordering Information

BAK-7195 BUTTONFIX[®] Fixation Button with Adjustable Loop