


Mini-Vim suture anchors

(Titanium)

The Mini-Vim suture anchors are titanium anchors available in 2.8mm and 3.5mm diameters and come preloaded with USP #2 BioFiber suture. In most cases, the implant can be manually inserted using the handled version. A drill is available for use in hard bone if desired. For arthroscopic applications, the Mini-Vim anchor can be inserted through the small diameter Mini-Vim spear eliminating the need for a cannula. This is ideal for SLAP and subscapularis repairs where a small stab incision and percutaneous delivery of the implant is preferred.

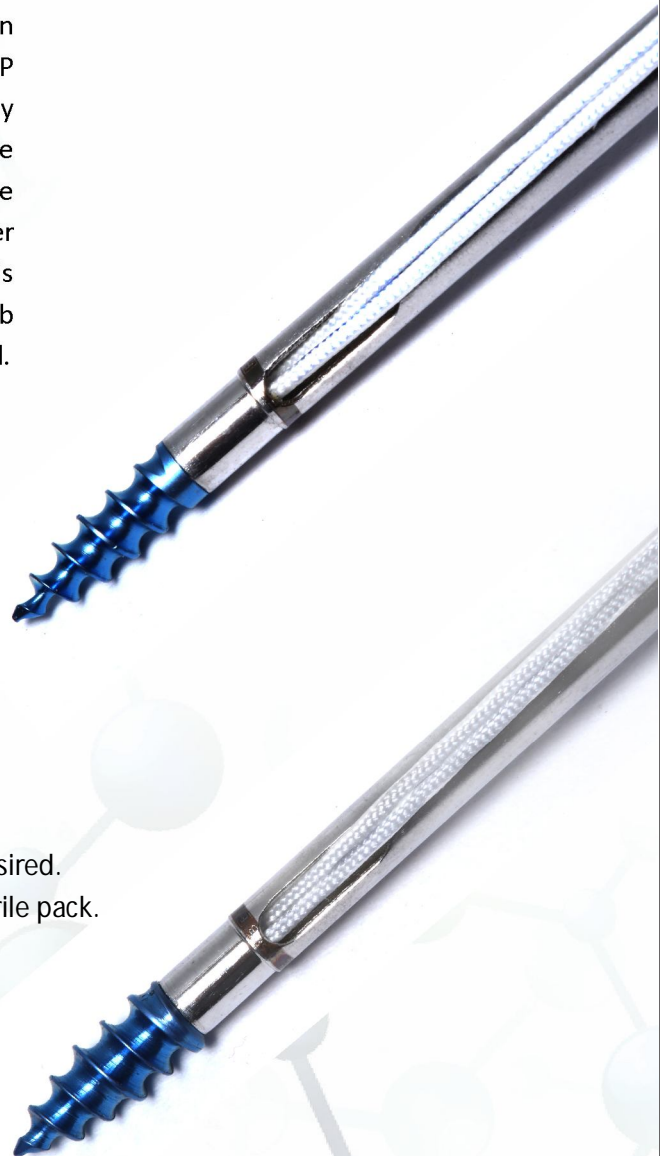

Key Features

- Tremendous pull-out strength.
- Self-drilling. A drill is available for use in hard bone if desired.
- Preloaded with One USP # 2 BioFiber in a single-use sterile pack.
- Also available in needled version.

Mini-Vim Suture Anchors

BAS-9001.28F	Mini-Vim suture anchor, dia.2.8mm,Ti (includes 1 no. BioFiber), Sterile
BAS-9040.28F	Mini-Vim suture anchor, dia.2.8mm, Ti (includes 1 no. BioFiber with needle), Sterile
BAS-9031.35F	Mini-Vim II suture anchor, dia.3.5mm,Ti (includes 1 no. BioFiber), Sterile
BAS-9083.35F	Mini-Vim II suture anchor, dia.3.5mm,Ti (includes 1 no. BioFiber with needle), Sterile

Related Instrumentation

BAS-9012	Mini-Vim Spear (Two-pronged) for BAS-9001.28F & BAS-9031.35F
BAS-9109	Mini-Vim Spear (Fishmouth) for BAS-9001.28F & BAS-9031.35F
BAS-9013	Trocar for BAS-9012 & BAS-9109
BAS-9026	Drill bit, dia.2.2mm for Mini-Vim (BAS-9001.28F)
BAS-9114	Drill bit, dia.2.8mm for Mini-Vim II (BAS-9031.35F)

BIOTEK - Chetan Meditech Pvt. Ltd.

Opp. V. S. Hospital, Ellisbridge, Ahmedabad-380 006. Gujarat, INDIA. Phone: +91 79 26578092
 Fax: +91 79 26577639 Email: info@biotekortho.com URL: www.biotekortho.com
 ISO 9001 • ISO 13485 • ISO 10002 • CE Certified • Indian FDA

Mini-Vim PK suture anchors (PEEK)

The Mini-Vim PK suture anchors are available in a variety of sizing and suture options. The Mini-Vim PK is a nonabsorbable PEEK suture anchor with a material eyelet which provides superior abrasion resistance due to PEEK's low coefficient of friction. Anchor insertion and delivery are made simple by drilling through a spear and inserting the anchor through the same spear into the drilled hole.

Mini-Vim PK Suture Anchors

BAS-9086.28F	Mini-Vim PK suture anchor, dia.2.8mm (loaded with one pc BioFiber), Sterile
BAS-9087.28F	Mini-Vim PK suture anchor, knotless, dia.2.8mm, Sterile
BAS-9088.30F	Mini-Vim PK suture anchor, dia.3.0mm (loaded with two pc BioFiber), Sterile
BAS-9089.35F	Mini-Vim II PK suture anchor, dia.3.5mm (loaded with two pc BioFiber), Sterile



Related Instrumentation

BAS-9110	Mini-Vim PK Spear (Two-pronged) for BAS-9088.30F & BAS-9089.35F
BAS-9097	Mini-Vim PK Spear (Fishmouth) for BAS-9088.30F & BAS-9089.35F
BAS-9098	Trocar for BAS-9097 & BAS-9110
BAS-9100	Drill bit, dia.2.4mm for Mini-Vim PK (BAS-9086.28F)
BAS-9101	Drill bit, dia.2.6mm for Mini-Vim PK (BAS-9088.30F)
BAS-9111	Drill bit, dia.2.6mm for Mini-Vim PK (BAS-9087.28F)
BAS-9102	Drill bit, dia.3.0mm for Mini-Vim II PK (BAS-9089.35F)
BAS-9012	Mini-Vim Spear (Two-pronged) for BAS-9086.28F & BAS-9087.28F
BAS-9109	Mini-Vim Spear (Fishmouth) for BAS-9086.28F & BAS-9087.28F
BAS-9013	Trocar for BAS-9012 & BAS-9109