

BIOTEK[®]
we put life in metal

BIOTEK[®]
we put life in metal



BIOTEK - Chetan Meditech Pvt. Ltd.


Plot No: MD 04, Near Teva Company, Sanand GIDC Gate No. 2,
Charal Gaam, Ahmedabad-382110. Gujarat, INDIA.
Email: contact@biotekortho.com Website: www.biotekortho.com

An ISO 13485 : 2016 Company
All Implants specified in the catalogue are CE certified



Influentia Designs-09374159337

Distributed by:

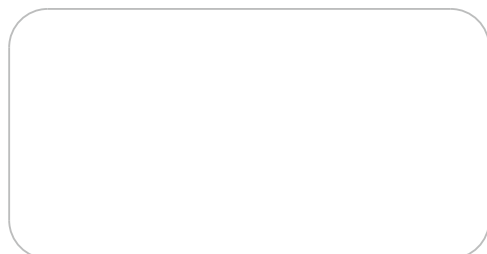

THIS MATERIAL DOES NOT CONSTITUTE MEDICAL ADVICE. IT IS INTENDED FOR INFORMATIONAL PURPOSES ONLY.
PLEASE CONSULT A PHYSICIAN FOR SPECIFIC TREATMENT RECOMMENDATIONS

BT/AR/11, Rev. 02

I n t r o d u c i n g

FIBERKNOT[®]
Soft Ligament Anchors

Small Hole...Secure Fixation...

www.biotekortho.com

FIBERKNOT® Ligament Anchors (Dia. 1.5mm and 1.8mm)

FIBERKNOT® Ligament Anchors (Dia. 2.9mm)

For Shoulder Instability Repair

FIBERKNOT® is an all suture anchor with smaller in size that does not mitigate strength. The FIBERKNOT® Soft Ligament Anchor is composed of ultra-high molecular weight polyethylene (UHMWPE) suture and is available in 1.5 mm (single loaded) and 1.8 mm (double loaded) sizes. FIBERKNOT® requires a smaller pilot hole compared to traditional anchors, minimizing bone removal and decreasing the overall anchor footprint, while still providing superior strength.

The FIBERKNOT® anchor is available with flexible driver for easy insertion through curved spear. It is available in various combination of fiber and fibertape.

The curved and straight spear of FIBERKNOT® enable delivery around a curve and perpendicular placement for posterior SLAP and inferior Bankart repairs.

The drill ensures accurate pilot hole placement with both straight and curved spears.



For Rotator Cuff Repair

The FIBERKNOT® helps overcome the traditional anchors for rotator cuff fixation methods with strong fixation, less bone removal and a simple, self-punching technique. FIBERKNOT® ligament anchor Dia. 2.9mm requires a pilot hole for the insertion. The punch is specially designed for FIBERKNOT® ligament anchor which provides the accurate hole placement for the anchor. The small 2.9mm footprint of FIBERKNOT® ligament anchor provides placement and fixation advantages when bone real estate is limited, especially during revision surgeries.

FIBERKNOT® is available in double loaded with various combination of fiber and fiber tape.




Ordering Information

Catalog No.	Product Description
10031	FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber, Dia. 2.9mm
10032	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber and 1 pc White-Blue Fiber Tape size 1.4, Dia. 2.9mm
10033	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber and 1 pc White-Black Fiber Tape size 1.4, Dia. 2.9mm
10034	FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber Tape size 1.4 (White-Blue & Solid Black), Dia. 2.9mm
10035	FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber Tape size 1.4 (White-Black & Solid Blue), Dia. 2.9mm
10052	FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber with needle, Dia. 2.9mm
10059	FIBERKNOT® Soft Ligament Anchor loaded with 3 pc Fiber, Dia. 2.9mm

Related Instrumentation

Catalog No.	Product Description
10055	Punch for FIBERKNOT® Soft Ligament Anchor, Dia. 2.9mm

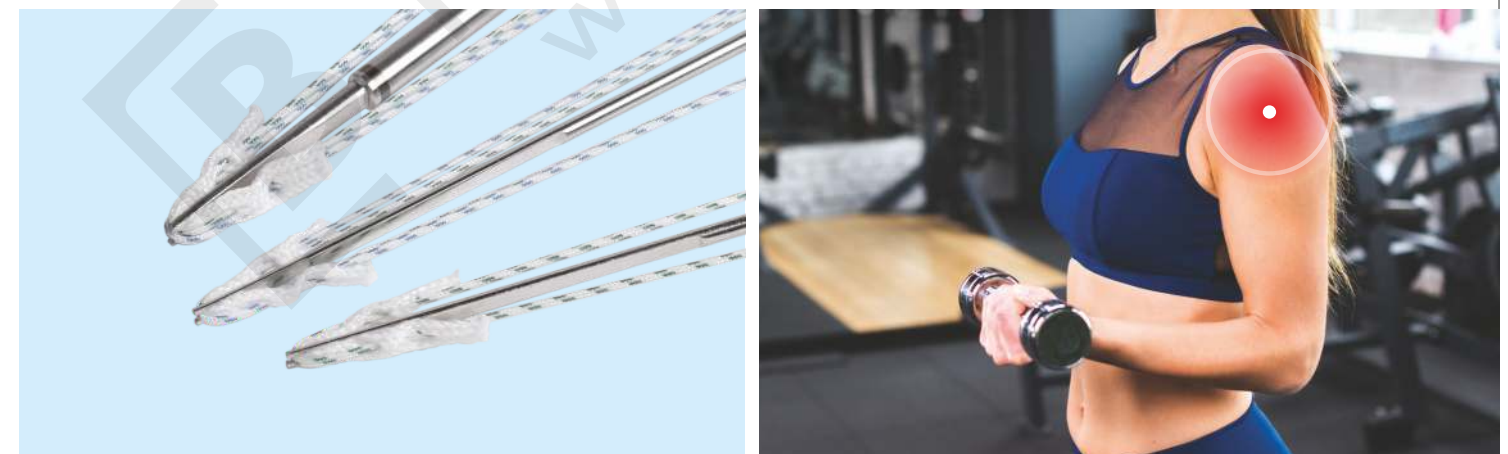

Ordering information

Catalog No.	Product Description
10025	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber with needle, Dia. 1.2mm
10026	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber, Dia. 1.5mm
10027	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber Tape size 1.4, Dia. 1.5mm
10028	FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber, Dia. 1.8mm
10029	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber and 1 pc Fiber Tape size 1.4, Dia. 1.8mm
10030	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc White-Black Fiber and 1 pc White-Blue Fiber Tape size 1.4, With Needle Dia. 1.8mm
10051	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber with needle, Dia. 1.5mm
10095	FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber with needle Dia. 1.8mm

Related Instrumentation

Catalog No.	Product Description
10040	Spear (Two-pronged), Straight for FIBERKNOT® Soft Ligament Anchor- Dia. 1.5mm & 1.8mm
10041	Spear (Two-pronged), Curved for FIBERKNOT® Soft Ligament Anchor- Dia. 1.5mm & 1.8mm
10042	Trocar, Straight for FIBERKNOT® Soft Ligament Anchor- Dia. 1.5mm & 1.8mm
10043	Trocar, Curved for FIBERKNOT® Soft Ligament Anchor- Dia. 1.5mm & 1.8mm
10045	Drill bit, Dia. 1.5mm for FIBERKNOT® Soft Ligament Anchor- Dia. 1.5mm & 1.8mm (SS)
10046	Drill bit, Dia. 1.5mm, for FIBERKNOT® Soft Ligament Anchor- Dia. 1.5mm & 1.8mm (Nitinol)