


Meniscal repair set

(Inside-Out Technique)

BIOCINCH[®]

BIOCINCH[®] Meniscal repair set comprises of two straight meniscal repair needles preloaded with USP size 2-0 or 0 suture. Meniscal repair needles are made of nitinol, and are ultrathin and strong, yet flexible for easy passage through guide sleeves. The suture is made of UHMWPE material provides ultra-high strength.

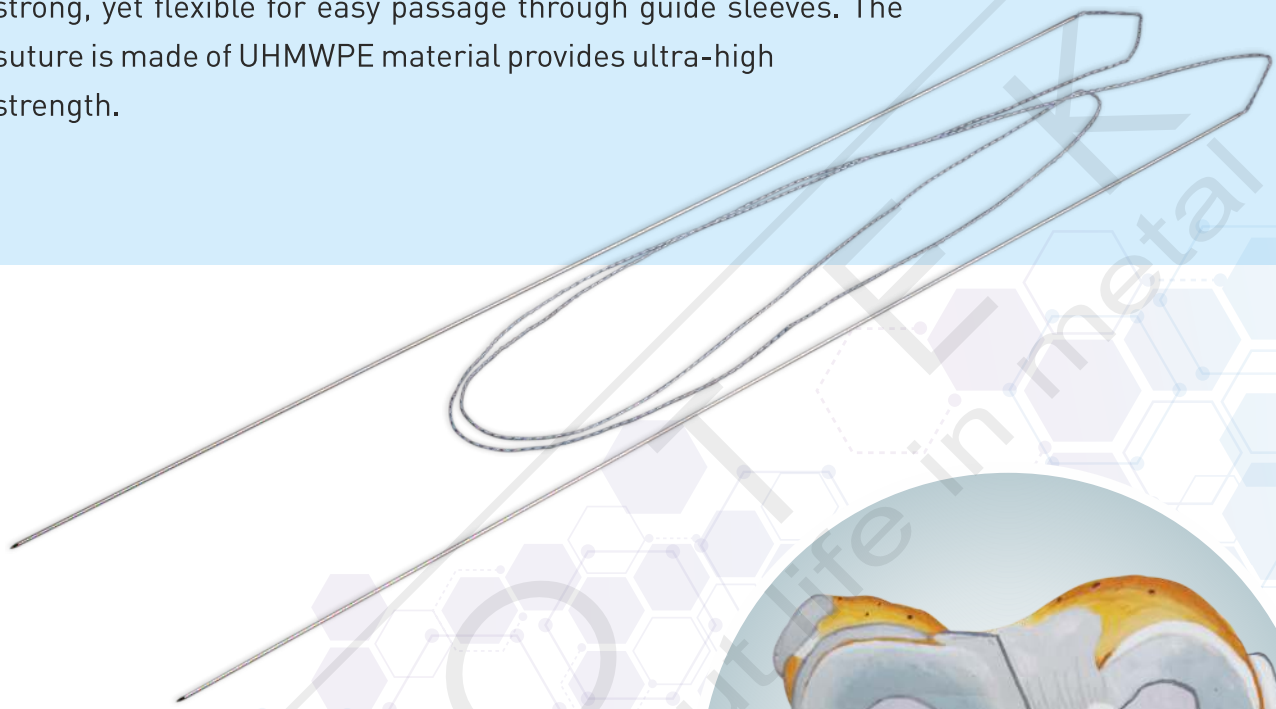


- Most frequently used technique.
- Simple to learn and replicate.
- Effective Suture placement-Vertical & Horizontal.

BIOCINCH[®] Meniscal repair set

BAK-7196 BIOCINCH[®] Meniscal repair set with needles; Fiber USP Size 2-0, suture Len. 75cm; straight needles-2 units, needle len. 25cm
BAK-7216 BIOCINCH[®] Meniscal repair set with needles; Fiber USP Size 0, suture Len. 75cm; straight needles-2 units, needle len. 25cm

BIOCINCH[®] Meniscal repair Instrument set

The BIOCINCH[®] Meniscal repair set has been designed to address meniscal tears in the anterior, middle and posterior zones of the meniscus using an "inside out" approach. The single lumen and double lumen pre-bent guide sleeves feature easy passage of needles, effective suture placement- vertical & horizontal, and to access easily on different zones of meniscus. The Curved Probe measures the width of the meniscus. The angled tip of the Meniscal Rasp is ideally shaped to access inside the meniscal tear for debridement prior to repair. Simple to learn and replicate.

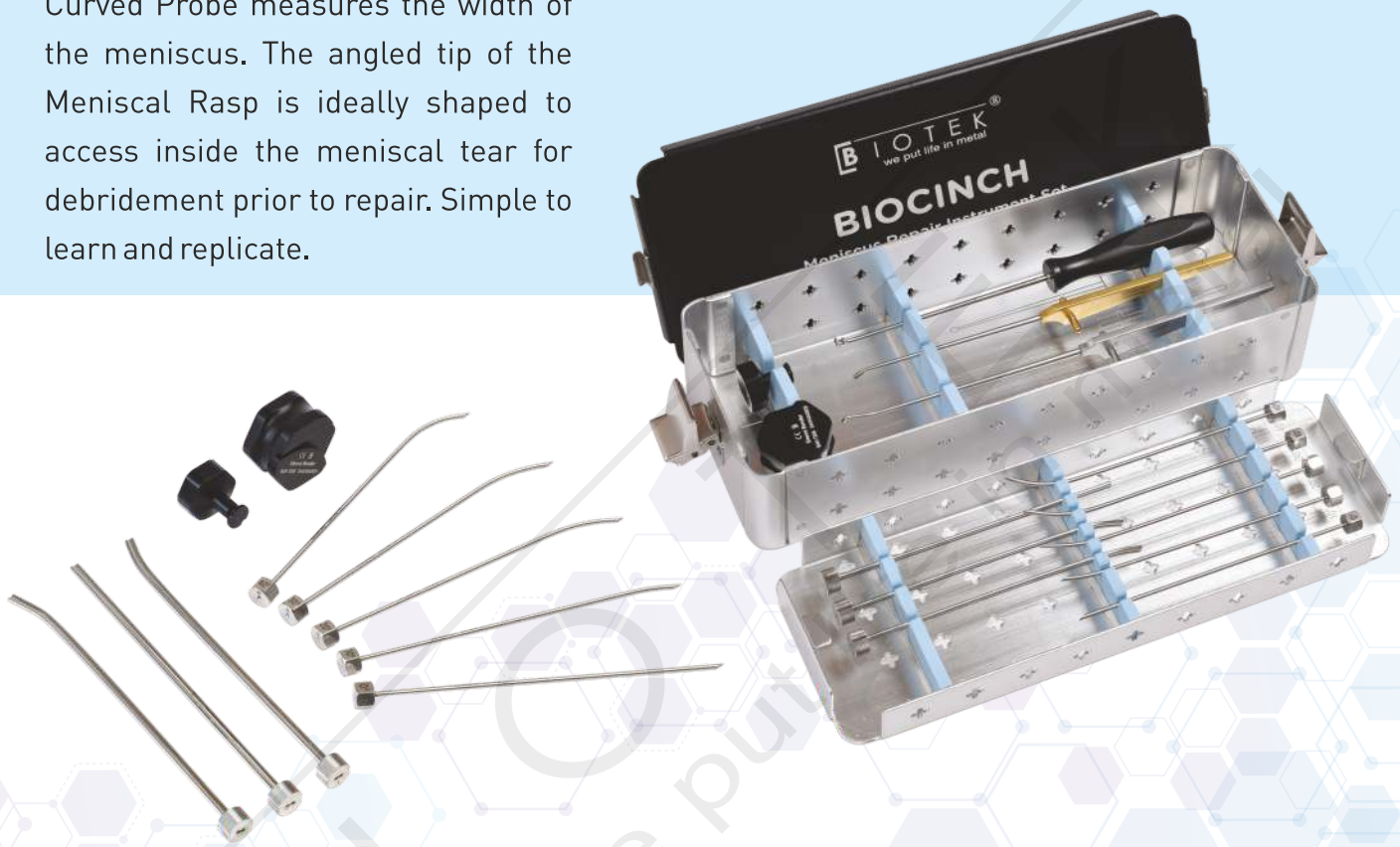

Related Instrumentation

BAK-7198	Curved Probe for Meniscus Repair
BAK-7199	Meniscus Rasp BAK-7204 Meniscus Currette
BAK-7200	Guide Sleeve (Double Lumen) for Meniscus Repair, Straight
BAK-7201	Guide Sleeve (Double Lumen) for Meniscus Repair, Curved, Left
BAK-7202	Guide Sleeve (Double Lumen) for Meniscus Repair, Curved, Right
BAK-7211	Guide Sleeve (Single Lumen) for Meniscus Repair, Straight
BAK-7212	Guide Sleeve (Single Lumen) for Meniscus Repair, Curved, 10°
BAK-7213	Guide Sleeve (Single Lumen) for Meniscus Repair, Curved, 20°
BAK-7214	Guide Sleeve (Single Lumen) for Meniscus Repair, Curved, 30°
BAK-7215	Guide Sleeve (Single Lumen) for Meniscus Repair, Curved, 40°
BAK-7205	Needle Pusher
BAK-7206	Sleeve Bender
BAK-7217	Graphic Case for Meniscus Repair Instruments

