

**BIOTEK**<sup>®</sup>  
we put life in metal

[www.biotekortho.com](http://www.biotekortho.com)



**ELBOW**  
SPORTS MEDICINE



Biotek is a global medical device company with a strong focus on new product development in the most advanced orthopedic sector: SPORTS MEDICINE

Biotek continues to experience unprecedented growth and demand for our products throughout the world; we currently have presence in all six continents around the world. However, Biotek remains a privately held company with a family business culture committed to delivering uncompromising quality to the health care professionals who use our products, and ultimately, the millions of patients whose lives we impact.

year founded  
**2003**



**40+**  
countries  
our products are  
used in



**100+**  
Sales partners  
globally



# MINI-VIM<sup>®</sup>

Ligament Anchor  
with Needle (Titanium)

ELBOW  
SPORTS MEDICINE



## MINI-VIM<sup>®</sup> Ligament Anchor with Needle (Titanium)

- The MINI-VIM<sup>®</sup> Ligament anchors are titanium anchors available in 2.8mm, 3.5mm diameter and come pre-loaded with USP #2 suture and needle.
- The implant can be manually inserted using the handled version. A drill is available for use in hard bone if desired.
- **Anchor Dimensions:** 2.8mm × 11.5mm, 3.5mm × 11.5mm  
**Type:** Screw in suture anchor  
**Suture Colour:** White

---

### ORDERING INFORMATION

---

BAS-9040.28F MINI-VIM<sup>®</sup> Ligament Anchor, Titanium loaded with 1 pc Fiber with needle, Dia. 2.8mm  
BAS-9083.35F MINI-VIM<sup>®</sup> Ligament Anchor, Titanium, loaded with 1 pc Fiber with needle, Dia. 3.5mm

---

### RELATED INSTRUMENTS

---

BAS-9026 Drill bit, Dia. 2.2mm for MINI-VIM<sup>®</sup> Ligament Anchor BAS-9001.28F & BAS-9040.28F  
BAS-9114 Drill bit, Dia. 2.8mm for MINI-VIM<sup>®</sup> Ligament Anchor BAS-9031.35F & BAS-9083.35F

# MINI-VIM PK®

Ligament Anchor  
(PEEK, Knotless)

ELBOW  
SPORTS MEDICINE



## MINI-VIM PK® Ligament Anchor (PEEK, Knotless)

- The MINI-VIM PK® Ligament anchor, Knotless are available in diameter of 2.8mm & 3.5mm.
- It is a non-absorbable PEEK suture anchor with a material eyelet which provides superior abrasion resistance due to PEEK's low coefficient friction.
- The knotless technology of MINI-VIM PK® Ligament anchor provides a step-saving alternative to conventional "knotted" suture anchors and eliminates "knot stacks" associated with soft tissue irritation.
- **Anchor Dimensions:** 2.8mm x 12.0mm, 3.5mm x 12.0mm  
**Type:** Tap-in/Bang-on suture anchor  
**Suture Colour:** White-Green (Holding Suture)

---

### ORDERING INFORMATION

---

BAS-9087.28F MINI-VIM PK® Ligament Anchor, PEEK, Knotless, Dia. 2.8mm  
BAS-9192.35F MINI-VIM PK® Ligament Anchor, PEEK, Knotless, Dia. 3.5mm

---

### RELATED INSTRUMENTS


---

BAS-9111 Drill bit, Dia.2.5mm for MINI-VIM PK® Ligament Anchor, Knotless-BAS-9087.28F  
BAS-9127 Sleeve for MINI-VIM PK® Ligament Anchor, Knotless-BAS-9087.28F

# FIBERKNOT®

Ligament Anchor  
with Needle (Soft)

ELBOW  
SPORTS MEDICINE

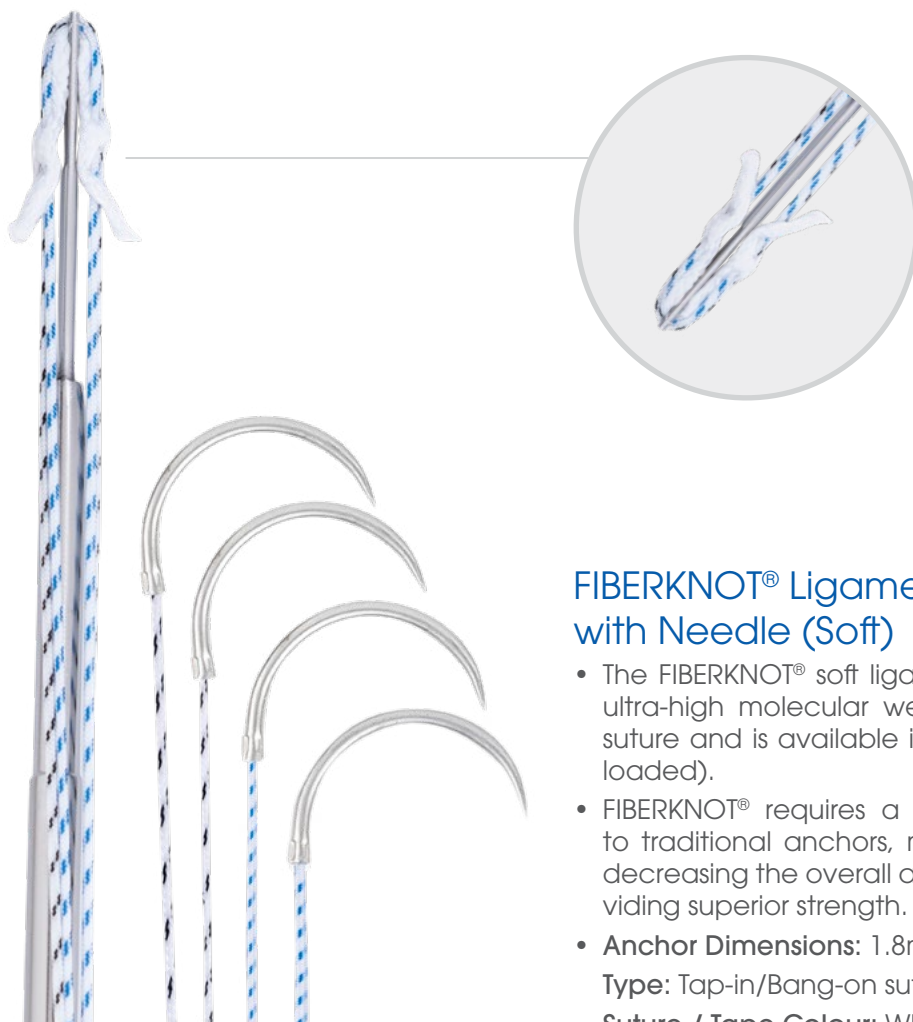

## FIBERKNOT® Ligament Anchors with Needle (Soft)

- The FIBERKNOT® soft ligament anchor is composed of ultra-high molecular weight polyethylene (UHMWPE) suture and is available in diameter of 1.8mm (double loaded).
- FIBERKNOT® requires a smaller pilot hole compared to traditional anchors, minimizing bone removal and decreasing the overall anchor footprint, while still providing superior strength.
- **Anchor Dimensions:** 1.8mm x 15.0mm  
**Type:** Tap-in/Bang-on suture anchor  
**Suture / Tape Colour:** White-Blue, White-Black

---

## ORDERING INFORMATION


---

- 10095 FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber with needle Dia. 1.8mm  
10030 FIBERKNOT® Soft Ligament Anchor loaded with 1 pc White-Black Fiber and 1 pc  
White-Blue Fiber Tape size 1.4, With Needle Dia. 1.8mm
- 

## RELATED INSTRUMENTS

---

- 10066 Drill Bit Dia 1.6mm , Len.- 100 mm for FIBERKNOT® Soft Ligament Anchor Dia 1.5mm and  
1.8mm with Needles
-

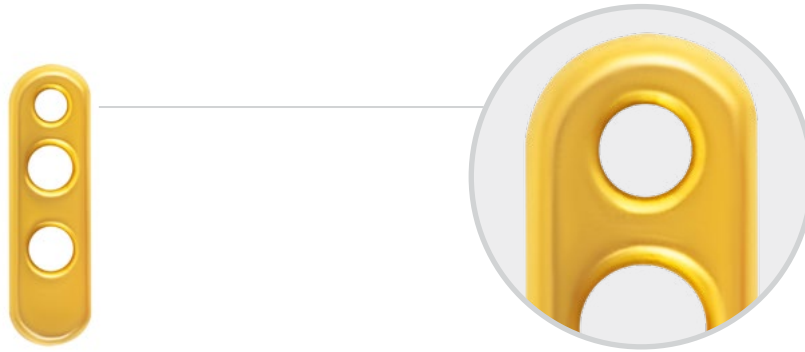

## AK-FIX Fixation Button

- The AK-Fix Fixation Button is a low profile Titanium button with slots enabling the insertion of the Biofiber loop.
- The oblong shaped button can be used with USP size 5 BioFiber Loop to treat Distal Biceps Repair.
- Button Dimensions: 10.0mm x 3.0mm x 1.3mm

---

### ORDERING INFORMATION

---

10104	AK-FIX Fixation Button
10122	AK-FIX Fixation Button, with Inserter

---

### RELATED INSTRUMENTS


---

10144	Suture Passing Inserter for AK-FIX Fixation Button
-------	--

# SOFTFIX PK®

Interference Screws,  
PEEK

## ELBOW SPORTS MEDICINE

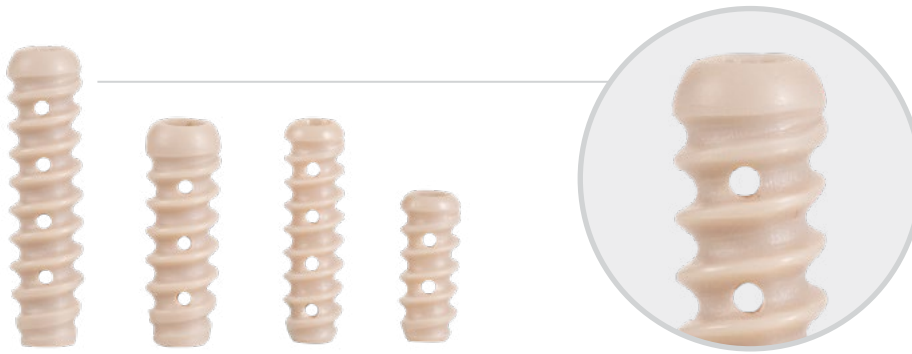

### SOFTFIX-PK® Interference Screws, PEEK

- Non-absorbable SOFTFIX-PK® interference screw incorporates PEEK-VESTAKEEP® polymer from Evonik® which exhibits an ideal combination of strength, stiffness and toughness that surpasses absorbable materials.
- It provides the interference screw fixation, while eliminating the need of trans osseous tunnels in tendon repairs and ligament reconstruction.
- It can be used for Distal Biceps Repair.
- Dimensions: Dia. 4.0mm x Len. 10.0mm & 15.0mm and Dia. 5.0mm x Len. 15.0mm & 20.0mm

---

#### ORDERING INFORMATION

---

10053	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 4.0mm, Len. 10mm
10054	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 4.0mm, Len. 15mm
BAK-7164.15	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 5.0mm, Len. 15mm
BAK-7164.20	SOFTFIX-PK® Interference screw, cannulated, roundhead; Dia. 5.0mm, Len. 20mm

---

#### RELATED INSTRUMENTS

---

10008	Screwdriver for SOFTFIX-PK® Interference Screw, Cannulated, Dia. 4.0mm
BAK-7223	Screwdriver for SOFTFIX-PK®, Dia. 5.0mm Interference Screws, Cannulated
10085	Fluted Reamer Cannulated, Dia. 4.0mm
10086	Fluted Reamer Cannulated, Dia. 5.0mm
BAL-024	Threaded Guidewire Dia. 1.2; Len. 180mm, SS
10094	Guide Wire, Blunt Tip for Dia. 4.0mm SOFTFIX-PK® Interference Screws Dia. 0.8mm; Len. 200mm (SS)
10096	Guide Wire, Blunt Tip for Dia. 5.0mm-8.0mm SOFTFIX-PK® Interference Screws Dia. 1.2mm; Len. 200mm (NITINOL)
10065	Drill Tip Suture Passing Pin, Dia. 2.4mm; Len. 200 mm



## BioFiber®

BioFiber® is a braided sterile surgical fiber prepared from fibres of ultra-high molecular weight long chain linear polyethylene. BioFiber® elicit a minimal of inflammatory reaction in tissues, followed by gradual encapsulation of the fiber by the fibrous tissue.

---

### ORDERING INFORMATION

BF-2W	BioFiber® Fiber USP Size 2 with MO-6 needle White, Len. 36"	BF-WWB	BioFiber® Fiber USP Size 2 Len. 36" Pack contains 2 units (White with needle & White-Blue with needle)
BF-2WC	BioFiber® Fiber USP Size 2 with MO-6 needle, White-Blue, Len. 36"	10068	BioFiber® Fiber USP Size 2-0 White, Len. 75cm
BF-2WK	BioFiber® Fiber USP Size 2 with MO-6 needle, White-Black, Len. 36"	10069	BioFiber® Fiber USP Size 0 White, Len. 75cm
BF-2WCK	BioFiber® Fiber USP Size 2 with MO-6 needle, White-Blue-Blue-Black, Len. 36"		

## BioFiber® Tape



## BioFiber® Tape

BioFiber® Tape is an ultra-high strength 2 mm width tape using a similar long chain polyethylene structure as the BioFiber® suture. The broad footprint of the BioFiber® Tape is appropriate for repairs in degenerative cuff tissue where tissue pull-through may be a concern.

---

### ORDERING INFORMATION

BF-2BST	BioFiber® Fiber Tape, size 2 White-Black, Len. 39"
BF-2ST	BioFiber® Fiber Tape, size 2, White-Blue, Len. 39"

## BioFiber® Loop



## BioFiber® Loop

BioFiber® Loop is a continuous loop of BioFiber® with straight needle. The needle is easy to handle and move freely on the suture to recenter itself after passing through tissue and facilitating even tension. Using BioFiber® Loop in conjunction with the Whipstitch technique drastically reduces the time spent and improves the fixation strength.

---

### ORDERING INFORMATION

BF-2WSB	BioFiber® Fiber Loop, Size 2 with straight needle, Len. 6cm White-Blue, Len. 24"
---------	--

---







**BIOTEK - CHETAN MEDITECH PVT. LTD.**  
MD 04, Sanand GIDC Gate No. 2,  
Ahmedabad - 382110. Gujarat, India.  
Email: [contact@biotekortho.com](mailto:contact@biotekortho.com)



DISTRIBUTED BY