

BIOTEK[®]
we put life in metal

www.biotekortho.com



HAND AND WRIST

Biotek is a global medical device company with a strong focus on new product development in the most advanced orthopedic sector: SPORTS MEDICINE

Biotek continues to experience unprecedented growth and demand for our products throughout the world; we currently have presence in all six continents around the world. However, Biotek remains a privately held company with a family business culture committed to delivering uncompromising quality to the health care professionals who use our products, and ultimately, the millions of patients whose lives we impact.

year founded
2003

40+
countries
our products are
used in

100+
Sales partners
globally

**SPORTS
MEDICINE**

MICRO-VIM™ Ligament Anchor (Titanium)

HAND AND WRIST



MICRO-VIM™ Ligament Anchor (Titanium)

- MICRO-VIM™ ligament anchor is made of titanium.
- The screw-in design of the anchor is ideal for soft-bone matrix.
- The MICRO-VIM™ Ligament Anchor is available in 1.8mm and 2.0mm diameters.
- The Inserter is preloaded with anchor and 1 pc USP size 2-0 suture with needle in dia. 1.8mm anchor and USP size 0 suture with needle in dia. 2.0mm anchor.
- MICRO-VIM™ ligament anchor can be used for carpal-metacarpal, tarsal-metatarsal and phalanges ligament repair of small bones.
- Anchor Dimensions : 1.8mm X 5.0mm, 2.0mm X 5.0mm
- Type: Screw in anchor
- Suture Colour : White-Blue
- Anchor Dimensions: Dia. 1.8mm x Len. 5.5mm; Dia. 2.0mm x Len. 5.5mm

ORDERING INFORMATION

BAS-9161.18F	MICRO-VIM™ Ligament Anchor, Titanium loaded with 1 pc Fiber USP 2-0 with needle, Dia. 1.8mm
BAS-9162.20F	MICRO-VIM™ Ligament Anchor, Titanium loaded with 1 pc Fiber USP 0 with needle, Dia. 2.0 mm

RELATED INSTRUMENTS

BAS-9150	Slotted Sleeve for MICRO-VIM™ & MICRO-VIM™ PK Ligament Anchor BAS-9161.18F, BAS-9162.20F, BAS-9163.18F & BAS-9164.20F
BAS-9149	Drill bit, Dia. 1.3mm for MICRO-VIM™ Ligament Anchor BAS-9161.18F
BAS-9148	Drill bit, Dia. 1.5mm for MICRO-VIM™ Ligament Anchor BAS-9162.20F

MICRO-VIM™ PK

Ligament Anchor
(PEEK)

HAND AND WRIST



MICRO-VIM™ PK Ligament Anchor (PEEK)

- MICRO-VIM™ PK ligament anchor is made of radiolucent PEEK polymer.
- The tap-in design of the anchor is ideal for hard-bone matrix.
- MICRO-VIM™ PK ligament anchor is available in 1.8mm and 2.0mm diameters. The Inserter is preloaded with anchor and 1 pc USP size 2-0 suture with needle in dia. 1.8mm anchor and USP size 0 suture with needle in dia. 2.0mm anchor
- MICRO-VIM™ PK ligament anchor can be used for carpal-metacarpal, tarsal-metatarsal and phalanges ligament repair of small bones.
- Anchor Dimensions: Dia. 1.8mm x Len. 8.0mm; Dia. 2.0mm x Len. 8.0mm
- Type: Tap in anchor
- Suture Colour: White-Blue

ORDERING INFORMATION

BAS-9163.18F MICRO-VIM™ PK Ligament Anchor, PEEK loaded with 1 pc Fiber USP 2-0 with needle, Dia. 1.8 mm
BAS-9164.20F MICRO-VIM™ PK Ligament Anchor, PEEK loaded with 1 pc Fiber USP 0 with needle, Dia. 2.0 mm


RELATED INSTRUMENTS

BAS-9146 Drill bit, Dia. 1.5mm for MICRO-VIM™ PK Ligament Anchor BAS-9163.18F
BAS-9146 Drill bit, Dia. 1.7mm for MICRO-VIM™ PK Ligament Anchor BAS-9164.20F
BAS-9150 Slotted Sleeve for MICRO-VIM™ & MICRO-VIM™ PK Ligament Anchor BAS-9161.18F, BAS-9162.20F, BAS-9163.18F & BAS-9164.20F

BUTTONFIX®-W

Two Fixation Buttons

HAND AND WRIST

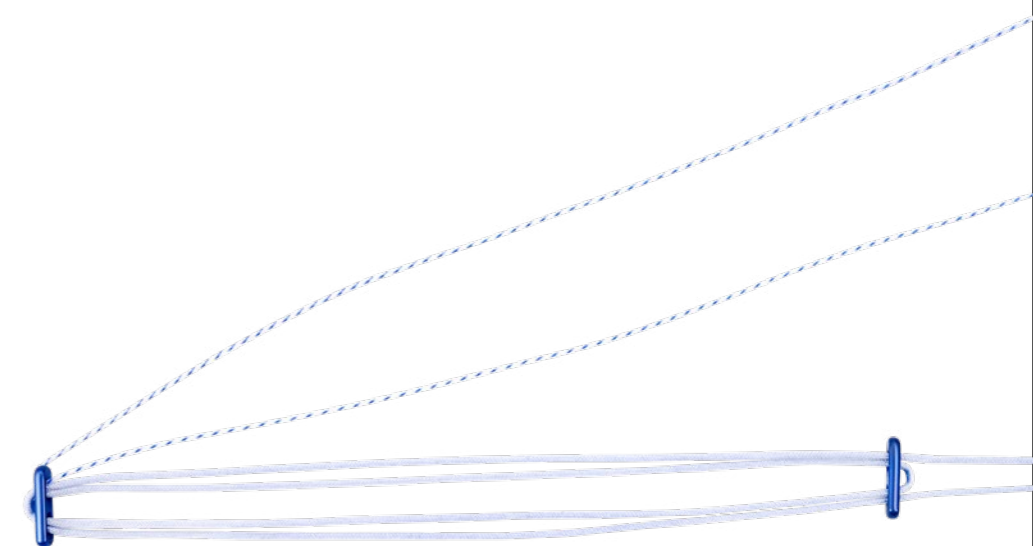

BUTTONFIX®-W (Two Fixation Buttons)

- BUTTONFIX®-W kit comprises of two oblong shaped fixation titanium buttons loaded with Fiber Loop USP size 2 and Fiber USP size 0 for flipping the button.
- This knotted adjustable loopsystem is ideal for the fixation of distal radioulnar joint (DRUJ) dislocation.
- This system is passed through the 2.7mm hole from the lateral to the medial side using a 2.4mm Suture Passing Pin with eyelet.
- This system comprises of a self-adjustable fiber loop size 2. The length of the loop can be adjusted by pulling the adjustable sutures one by one and final knotting required..
- Button Dimensions: 2.5mm x 1.1mm x 8.0mm

ORDERING INFORMATION

10081 BUTTONFIX®-W Two Fixation Buttons loaded with Fiber Loop USP 2 and Flipping Fiber USP 0

RELATED INSTRUMENTS

BAL-024 Threaded Guidewire Dia. 1.2; Len. 180mm, SS
BAL-014 Drill Bit, Dia: 2.7mm, Cannulated
10065 Drill Tip Suture Passing Pin, Dia. 2.4mm; Len. 200 mm
10010 Suture Passing wire, Blunt Tip with Loop Len. 200mm

BUTTONFIX®-W

One Fixation Button

HAND AND WRIST



BUTTONFIX®-W (One Fixation Button)

- BUTTONFIX®-W kit comprises of one oblong shaped fixation titanium button, preloaded with Fiber Loop USP size 2.
- The system requires a very small hole of 1.5mm to pass from radial to ulnar side metacarpal through the 1.5mm tapered trocar tip Suture Passing Wire with loop.
- This system supports and maintains the thumb and index metacarpals in proper alignment after trapezial resection or removal for treatment of CMC arthritis and thumb CMC instability.
- WT-Fix fixation button is used with this system for final fixation on the ulnar side of the index finger after the preloaded button is pulled from the radial side of the first metacarpal.
- This knotted system is ideal for the treatment of CMC Arthritis and/or CMC instability.
- Button Dimensions: 2.5mm x 1.1mm x 8.0mm

ORDERING INFORMATION

10204	BUTTONFIX®-W One Fixation Button loaded with Fiber Loop USP 2
10105	WT-FIX Fixation Button

RELATED INSTRUMENTS

10224	Suture Passing wire, Trocar Tip with Loop Dia. 1.5mm, Len. 250mm
-------	--

FIBERKNOT®

Soft Ligament Anchors

HAND AND WRIST



FIBERKNOT® Soft Ligament Anchors

- FIBERKNOT® is an all suture anchor and is composed of ultra-high molecular weight polyethylene (UHMWPE) suture.
- It is available in 1.5 mm (single loaded) and 1.8 mm (double loaded) sizes with needle.
- FIBERKNOT® requires a smaller pilot hole compared to traditional anchors, minimizing bone removal and decreasing the overall anchor footprint, while still providing superior strength.
- The FIBERKNOT® anchors can be used for any soft tissue fixation, Eg.: CMC Arthritis, TFCC Repair, SL Repair, etc.
- Small inserter provides easy and quick insertion.
- Anchor Dimensions: Dia. 1.5mm x Len. 15.0mm; Dia. 1.8mm x Len. 15mm
- Type: Tap in anchor

ORDERING INFORMATION

10051	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc Fiber with needle, Dia. 1.5mm
10095	FIBERKNOT® Soft Ligament Anchor loaded with 2 pc Fiber with needle Dia. 1.8mm
10030	FIBERKNOT® Soft Ligament Anchor loaded with 1 pc White-Black Fiber and 1 pc White-Blue Fiber Tape size 1.4, With Needle Dia. 1.8mm

RELATED INSTRUMENTS

10006	Slotted Sleeve for Foot & Ankle Drill bits-10004, 10005, 10067
10066	Drill Bit Dia 1.6mm, Len.- 100 mm for FIBERKNOT® Soft Ligament Anchor Dia 1.5mm and 1.8mm with Needles



BIOTEK - CHETAN MEDITECH PVT. LTD.

MD-04, Sanand GIDC-II,
Ahmedabad, Gujarat-382110. India.
Email: contact@biotekortho.com



DISTRIBUTED BY